

# OLYMPUS

## TD-20 MANUAL DISINFECTOR



A further product in a design  
co-ordinated endoscopy support range




# TD-20 MANUAL DISINFECTOR

The TD-20 Manual Disinfector, allows the cleaning and disinfection of Olympus flexible endoscopes.

## FEATURES INCLUDE:

- **Framework design** – co-ordinated with the popular Olympus workstation range.
- **Integral hand pump** – situated in the trolley leg for improved operation and aesthetics.
- **Three fluid containers with lids** – to allow total scope processing from cleaning to disinfection and final rinsing. The specially designed container lids minimise fluid evaporation, while the integral pouring lip and hand grip on the containers facilitate filling with little or no splashing.

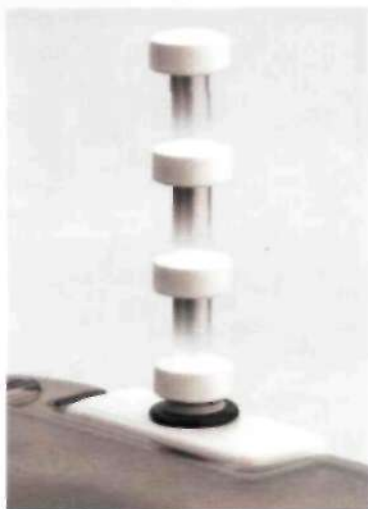

- **Fluid irrigation connectors** – all connectors are integral to the top tray, minimising the risk of splashing.
- **Smoked polycarbonate lid** – helps prevent disinfectant evaporation while allowing visual monitoring of the scope at any stage of the cycle.
- **Integral timer with alarm** – allows the disinfectant contact time to be accurately monitored.
- **Instrument tray** – purpose designed to accept Olympus OES fiberscopes and EVIS videoscopes in the minimum amount of disinfectant.

## CONSTRUCTION:

The TD-20 Manual Disinfector has a powder-coated steel frame, vacuum-moulded ABS plastic instrument tray and shelf and smoked polycarbonate lid. Three buckets with lids are supplied, together with channel connectors for upper and lower GI OES fiberscopes and EVIS videoscopes.

## OCCUPATIONAL HEALTH:

The Manual Disinfector requires the use of a cold fluid chemical disinfectant. It is important that whichever disinfectant is selected, the appropriate health and safety requirements for ventilation and personal protection are observed. Where appropriate the requirements of the COSHH regulations (UK) must be complied with.

## SPECIFICATION:

Dimensions:	960(W) x 930(H) x 510(D) mm
Weight:	27 kg
Instrument Tray Capacity:	6 litres



Olympus reserves the right to alter the specification without notice.

## OLYMPUS®

OLYMPUS OPTICAL CO., LTD.

San-Ei Building, 22-2, Nishi-Shinjuku 1-chome, Shinjuku-ku, Tokyo, Japan

OLYMPUS KEYMED

KeyMed House, Stock Road, Southend-on-Sea, Essex S52 5QH, United Kingdom

OLYMPUS OPTICAL CO., (EUROPA) GMBH.

(Premises-Goods delivery) Wendenstraße 14-16, 20097 Hamburg, Germany  
(Letters) Postfach 10 49 08, D-20034 Hamburg, Germany

OLYMPUS AMERICA INC.

Two Corporate Center Drive, Melville, New York 11747-3157, USA

OLYMPUS SINGAPORE PTE LTD.

401B River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373, Singapore



ISO 9001 Registered  
Certificate No. FM 20993